8. Alarm Tone(s)

▲Press button **3** first, then press button **1** to chose the next alarm tone.

9. Testing

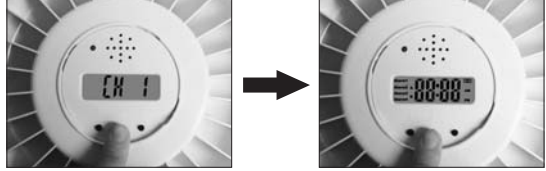
▲press button **1** for about 12 seconds, it display CH

The Automatic Pill Dispenser has 3 different alarm tones and also silence to accommodate different hearing levels. Symbols A, B and A+B on the display indicate the selected alarm tone, if "silence" is selected, then neither A nor B appears on the display. To change the alarm tone, first hold down button 3 and then, at the same time press and hold down button 1 to choose the next alarm tone. Continue to hold button 3 and press button 1 to listen to and select the different tones.

CH Test Mode -When the display shows CH, pill box is ready to be tested. To perform a system function test, press button 2. The unit will precede through the following system test modes.

■ CH 1 Display Test

The display shows CH 1. All segments on the display shall light up.

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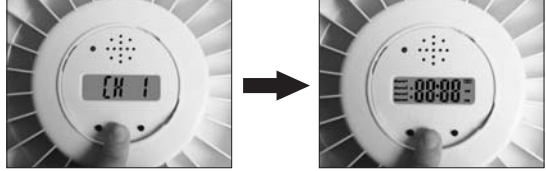
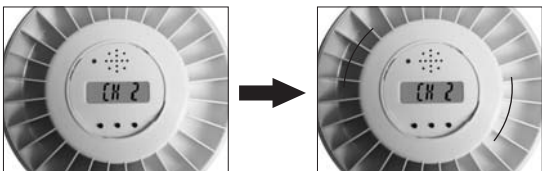
■ CH 2 Display Test

The display shows CH 2. The motor shall rotate the tray one compartment.

■ CH 3 Alarm Signal Test

The display shows CH 3. The alarm should sound and LED light lights up. Turn the pill box over to test RESET function (as to dispense medicine). The sound and LED light should be stopped. The CH 4 message will appear on the display when the pill box is not turned over within 10 seconds.

The light becomes red

●If the Automatic Pill Dispenser does not rotate properly the motor will stop after a few seconds and display the error message "CH6". Exit the test mode by pressing button 3.

●If there is a malfunction, please replace the batteries to reset.

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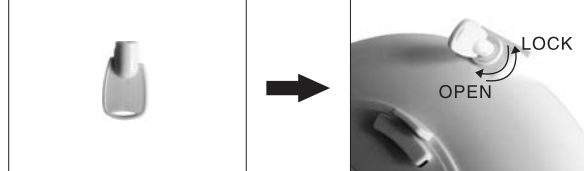
10. How to Lock and Unlock the Automatic Pill Dispenser

Insert key and turn to left 1/4 turn to lock, Insert key and turn to right 1/4 turn to unlock.

11. How to Stop Alarm

To stop the alarm and dispense the pills, pick up the Automatic Pill Dispenser and tilt over, dumping the pills into the hand or a suitable container.

12.The Low Battery Alert

When the Low Battery indicator (as picture) appears on the LCD, the batteries are low and need replacing.

Distributed by
C&A Marketing, Inc.
Ridgefield Park, NJ 07660
info@myivation.com
1-866-949-3049

MADE IN CHINA
IVATION is a trademark of
C&A IP Holdings, LLC
©2015 All Rights Reserved

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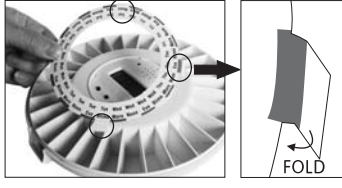
The Automatic Pill Dispenser is supplied with 4 discs. The discs identify the day of the week, and the appropriate number of doses per day (1, 2, 3, or 4 doses per day).

The appropriate disc for the medication regime is then placed on the tray; fold tabs down to secure in pill compartment, prior to filling it, and used as a guide as to where to place the pills.

Organizing the medication, assigning dose times and determining what medications go into each slot is the responsibility of you, your caregiver, and/or your health care provider.

The refill frequency will depend upon the number of daily medications for example:
One time a day medication refill every 28 days
Twice a day medication refill every 14 days
Three times a day medication refill every 9 days
Four times a day medication refill every 7 days

5. 12/24-Hour Time Formats

▲ Press button [2] and button [3] simultaneously

The automatic pill dispenser automatically defaults to a 24 hour time format. To switch to a 12 hour format, press buttons 2 and 3 simultaneously. The AM or PM symbol will appear on the right side of the display when the clock is in the 12 hour time format. Release the buttons when the desired time format is displayed.

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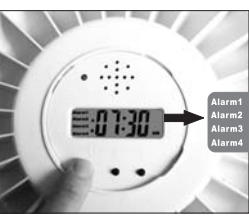
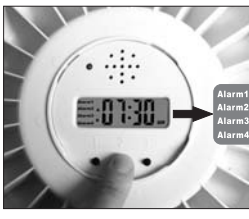
(button 1- left; button 2- middle; button 3- right)

▲ Press button [1]

1. Press and hold button 1 for a few seconds (until the symbol ":" blinks) and then release.

BLINK

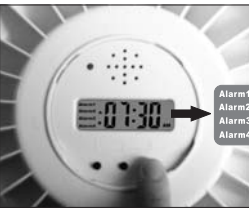
▲ Press button [2]

2. Set the correct time by pressing button 2 (hours) and button 3 (minutes).

BLINK

Set hour

▲ Press button [3]

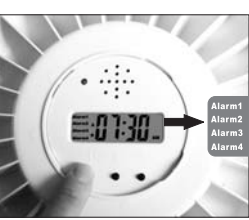
3. Store clock time by pressing button 1.

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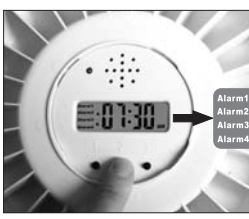
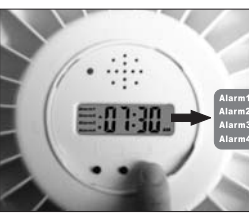
CLIGNOTEMENT

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7. Setting and Checking the Alarm Times

1. Press button 1 and then release. The alarm 1 symbol will begin to blink. To set alarm time 1, press buttons 2 (hours) and 3 (minutes). If you wish to cancel the alarm, press button 2 until is displayed, which occurs after 23:00 or 11:00 p.m.

2. After you have programmed the correct alarm time, press button 1 again to store the alarm time and to advance to the next Alarm setting.

3. The Alarm 2 symbol will blink. To set alarm time 2, use buttons 2 & 3.

4. Continue as above to set Alarm setting 3 and Alarm setting 4.

5. End by pressing button 1 until the Alarm symbols cease to blink.

■ If buttons 2 or 3 are depressed, time display progresses automatically. If you let the display continue to blink in the setting mode, it will automatically return to normal function after one minute.

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