VG832A
Deluxe Overbed Table

PLEASE NOTE:

TOOLS FOR ASSEMBLY ARE INCLUDED BUT FOR EASIER ASSEMBLY OTHER TOOLS MAY BE REQUIRED.

SOME SKILL AT ASSEMBLING REQUIRED

This file is available to view and download as a PDF at www.aidapt.co.uk. Sight impaired customers can use a free PDF Reader (such as adobe.com/reader) to zoom in and increase the text size for improved readability.

Fixing and Maintenance Instructions

NB. This equipment must be installed by a competent person.

VG832A Self Assembly with castors

Weight Limit: 20kg

Do not exceed the weight limit on the table as this could damage the equipment or put the user at risk. Exceeding the weight limit may also cause the table to tilt.

Do not lean on the table or use it as a support aid when rising from a chair or bed.

If the table is to be used with hot items (eg hot drink) then a risk assessment is advised as to scalding etc. because part of the table can tilt.


BEFORE USE

Carefully remove all packaging and avoid using knives or other sharp instruments as this may damage the surface of the equipment. Check the product for any visible damage and if you see any damage or suspect a fault do not use the equipment and contact your supplier.

ASSEMBLY

1. Attach top of Column (photo C) to the 'plate'on the side table top (photo E) using 4 x nut, bolt and washer. The nuts are 'nylocs' and are a safety device to stop the nut and bolt on screwing. The nuts and bolts need to be fully tightened using the spanner and allen key and the nyloc will give resistance all the way until it is fully tightened. The nuts are on the inside of the plate with the bolt head showing on the outside of the plate.

2. Attach the Base (photo B) to the bottom of the Column (photo C) using 1 x bolt and 1x nut, bolt and washer (again nyloc) making sure completely tightened.

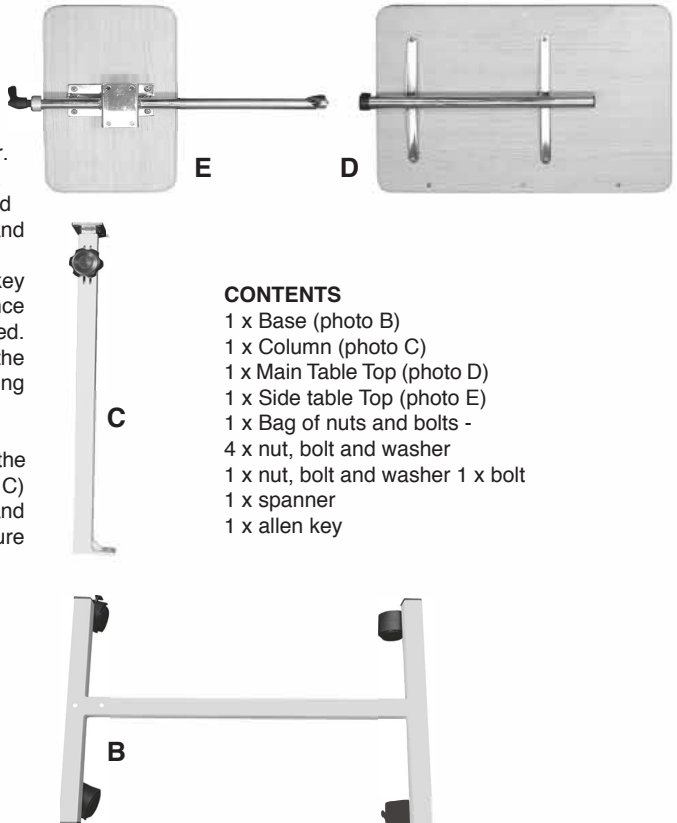
CONTENTS

- 1 x Base (photo B)
- 1 x Column (photo C)
- 1 x Main Table Top (photo D)
- 1 x Side table Top (photo E)
- 1 x Bag of nuts and bolts -
- 4 x nut, bolt and washer
- 1 x nut, bolt and washer 1 x bolt
- 1 x spanner
- 1 x allen key

ASSEMBLY (CONTINUED)

3. The Main Table Top (photo D) is then fitted on to the Side Table Top (photo E). When fitted it as par Photo F.

4. Make sure the tightening lever (photo G) and height adjustment knob (photo H) are tight so that the Main Table Top will not tilt and the table top will not slip downwards

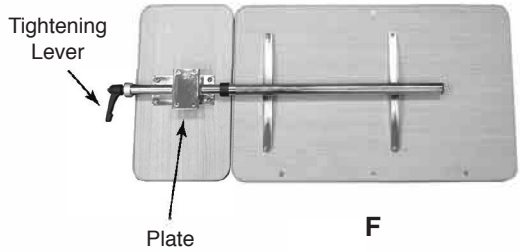USAGE

The table can be used in several ways.


First of all it can be used as Photo A. The Main Table Top can then be tilted as per Photo G - this is done by undoing the tightening nut, then tilting the Main Table Top to the required tilt and then retightening the nut.

The Main Table Top can also be extended sideways again by loosening the tightening nut, then extending the table and re tightening the nut (photo H). If the table is extended please make sure it is safe to do so in the particular circumstances (risk assessment) but would recommend that the table width of the whole unit should not exceed 1 metre.


NB. Please make sure the tightening nut is tight so that the Main Table Top will not tilt from its set position. The Side Table Top will not tilt and is therefore better suited for hot drinks etc so as to avoid spilling/scalding.

HEIGHT ADJUSTMENT

The height of the table can be adjusted by first loosening the height adjustment knob and then selecting the required height. Re tighten the height adjustment knob making sure it is tight enough so that the table will not slip downwards. It is recommended that the maximum height of the table is 1 metre.

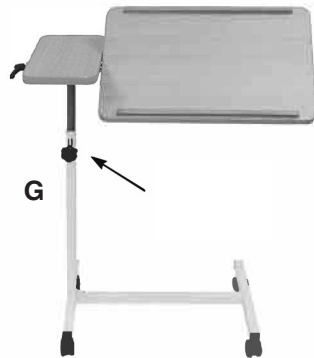DIMENSIONS

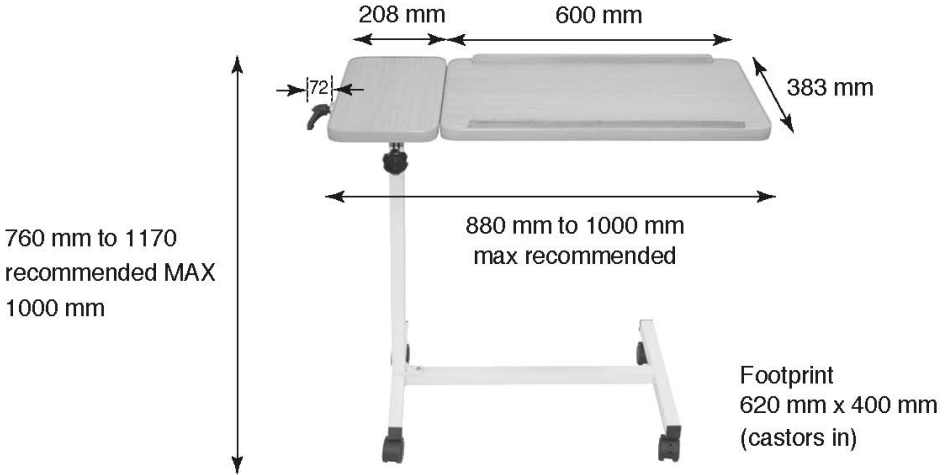CLEANING

A useful guide to decontamination is MHRA leaflet DB2003 (06)

Do not use abrasive cleaners as this could damage the surface of the equipment. For cleaning, wipe with a clean cloth soaked in detergent solution and wring out. Depending on the detergent manufacturer's instructions the solution may be heated.

Wipe all surfaces of the equipment and the cloth should be re-immersed periodically and the detergent changed frequently. The detergent should be rinsed or wiped off using a clean cloth immersed in clean water and wrung out.

To disinfect make sure the equipment has been cleaned and dried and then wipe over with a clean cloth and a compatible liquid disinfectant. After disinfecting clean with clean water and dry.

When using detergent or disinfectants then please refer to manufacturers' instructions and if unsure then consult your local Infection Control Specialist for advice.

The metal frame of this equipment can be cleaned/disinfected by moist heat at 80oC for a minimum of 1 minute or autoclaved at 80oC for at least 1 minute. After cleaning or disinfecting treat all joints with a water dispersing agent, e.g. WD40.

NB. The wooden table top is not suitable for cleaning with moist heat or autoclaving.

REISSUE

If you reissue or are about to reissue this product, please thoroughly check all components for their safety. If in any doubt, please do not issue or use, but immediately contact your supplier for service support.

CARE & MAINTENANCE

Please make a safety check of the product at regular intervals or if you have any concerns.

SPECIFICATIONS

Net Weight.....	9.9 kg
Gross Weight.....	10.8 kg
Dimensions of outer packing box	885 x 140 x 420 mm

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Aidapt Bathrooms Ltd.