Version Mars 2026

IDENTITÉS

CE

FAUTEUIL RELEVEUR MYST

LIFT CHAIR MYST

REF :
821160.TISSU
821160.CUIR

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INTRODUCTION :

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
INSTRUCTION MANUAL
ENG

Safety precautions:

Maximum user weight: 130 kg. Do not use the chair if the user's weight exceeds this limit in order to avoid any risk of malfunction or injury.

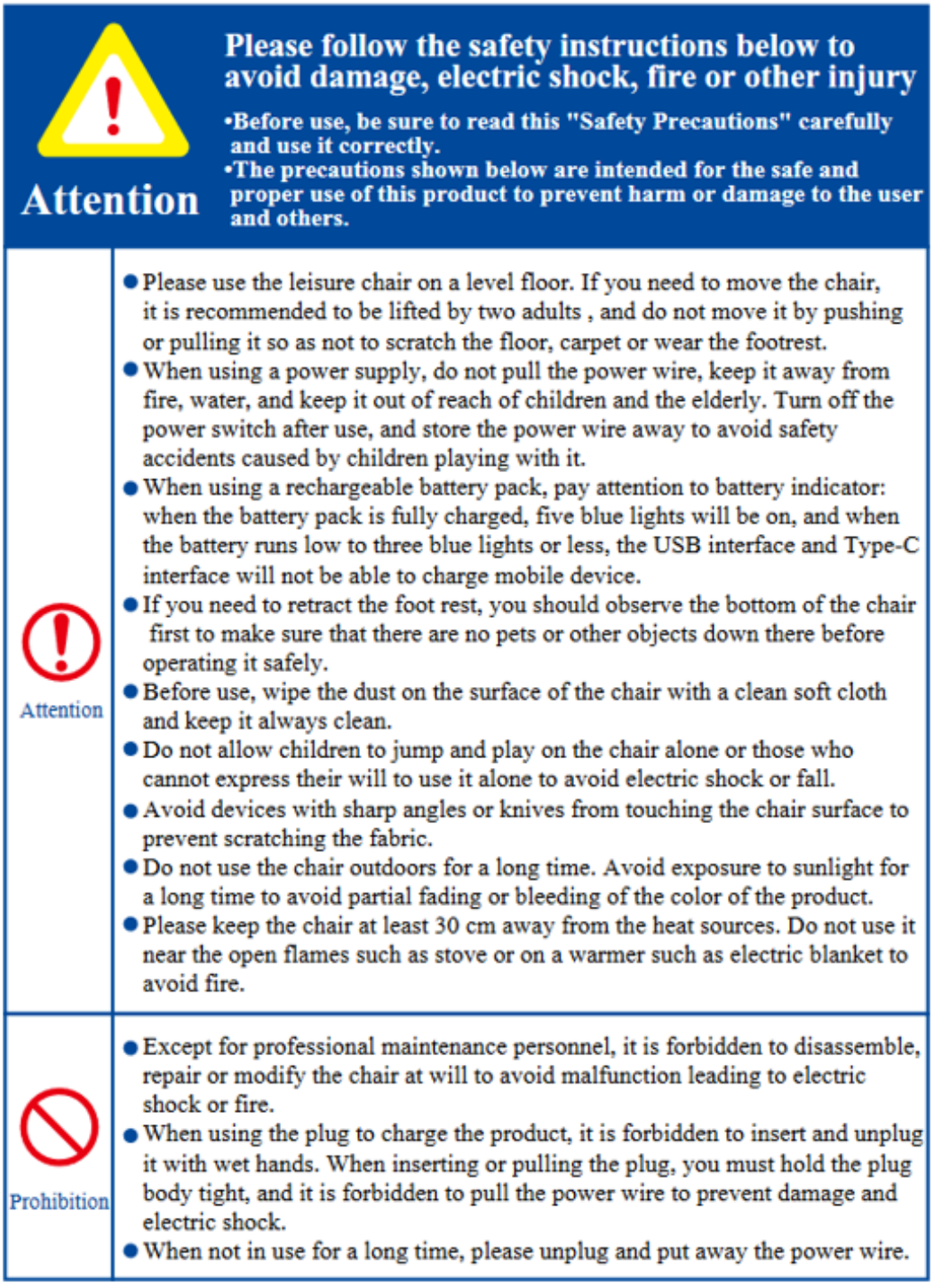Installation Method:

Components include:

❶ Back cushion X1 ☑


❷ Seat X1 ☑

❸ Power supply X1 ☑


Installation(Assembly) procedure description:

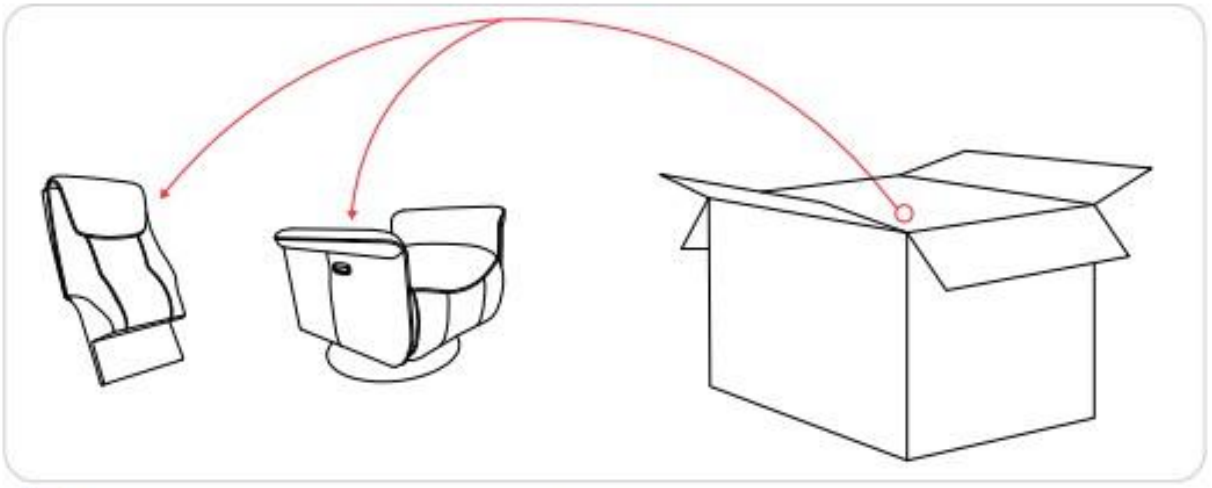1 Remove the back cushion and seat from the carton, lay it flat on the ground。

2 Align the locks of back cushion with the latch of the seat and press the latch firmly. When you hear a"click", the installation is successful.

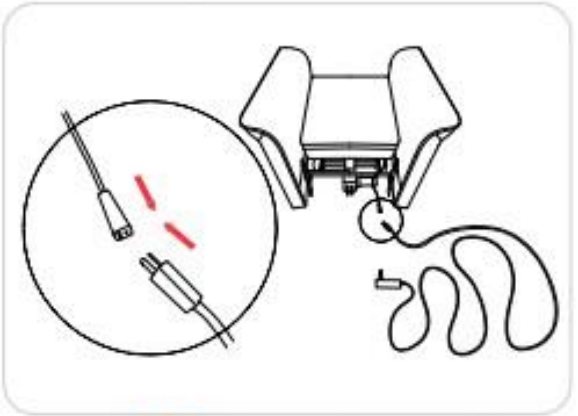3 3.Connect the control box to the external power cable.

Environment in use and Application Method:


Environment in Use:

It is recommended to place it in a flat place, with a space around it that can allow the chair back and foot rest to unfold, so as to avoid collision.

When lowering the back rest or extending foot rest, in order to avoid collision with surrounding objects, a space of more than 45 cm in front and 42 cm in the rear should be reserved in advance.

Do not use the product outdoors. Avoid exposure to sunlight for a long time to avoid local fading or blooding of colors of the product and rapid aging of the fabric.

Please keep the chair at least 30 cm away from the heat sources. Do not use it near the open flames such as stove or on a warmer such as electric blanket to avoid fire.

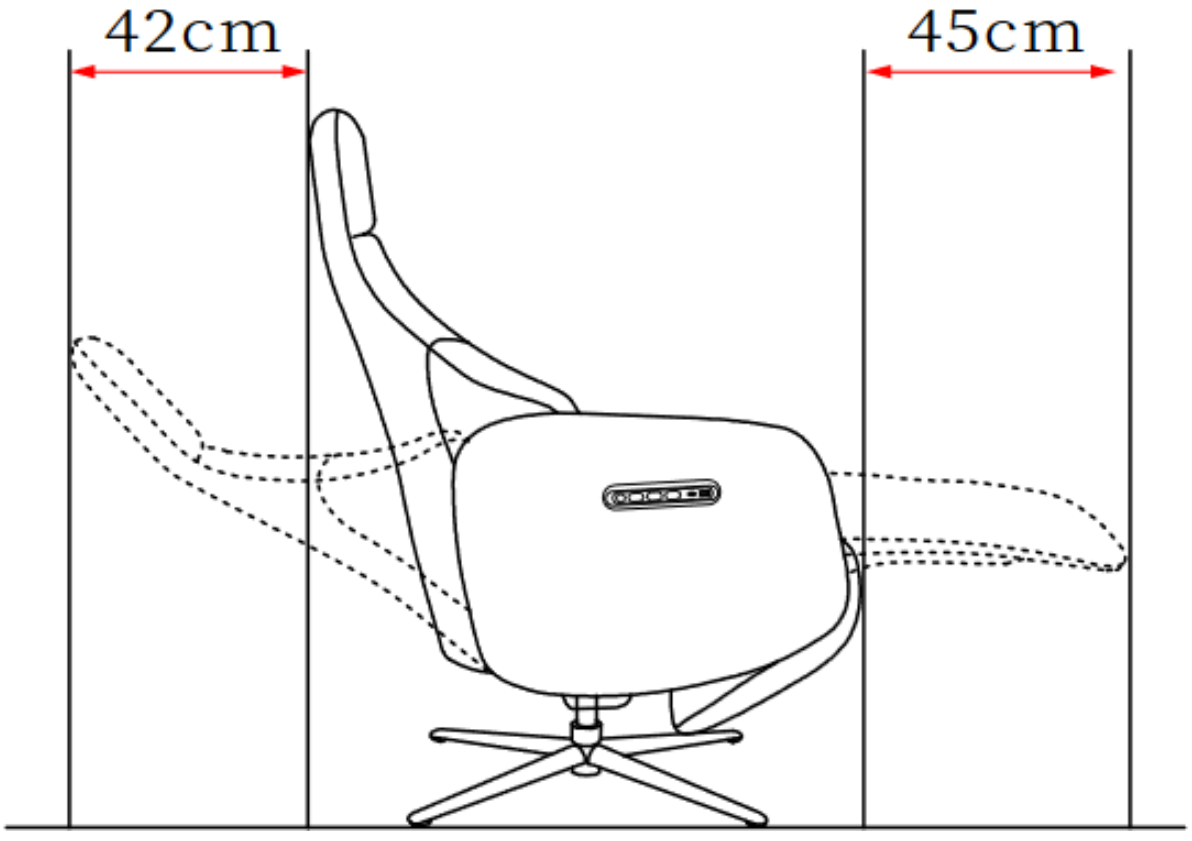
Environment in use and Application Method:


Application Method:

·Use method of the switch button [applicable to the Manual headrest function]

The switch buttons ①/② are used to control the opening/retraction of the backrest, the buttons ③/④ are used to control the footrest opening/retraction+seat lift .
⑤ is USB interface

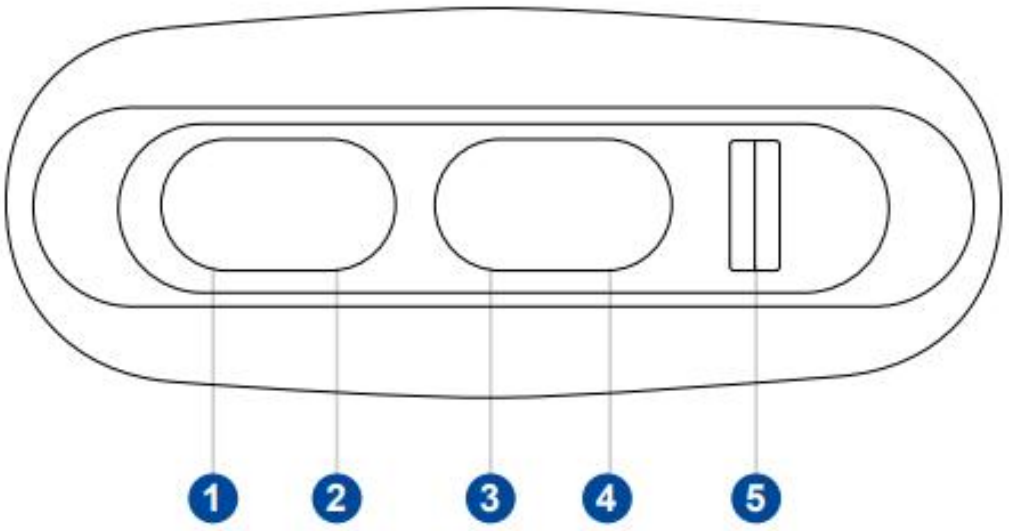1. Backrest leaning
2. Backrest retraction
3. Footrest retraction+seat lift
4. Footrest opening
5. USB interface

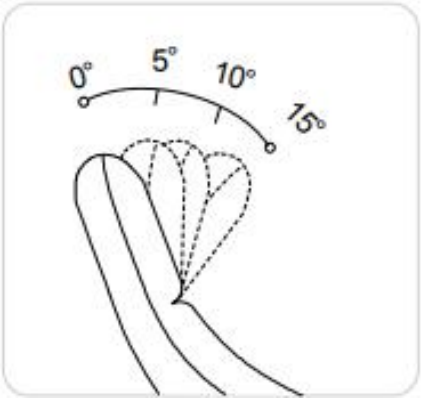·Manual headrest angle adjustment

The headrest mechanism of this product has 4 levels. When operating, grasp both sides of the headrest with both hands and push forward with a little force. When you hear a "click" or feel slight resistance, you can stop at the current position: let go of it when the maximum angle is reached, then it can be automatically restored to the initial angle.


Attention:


- Clean the chair surface with a clean soft cloth before use, and keep it clean frequently.
- Teach children not to jump on the seat surface of the chair, so as not to deform after partial compression and affect the use.
- Avoid devices with sharp angles or knives from touching the chair surface to prevent scratching the fabric.
- If you need to retract the foot rest, you should observe the bottom of the chair first to make sure that there are no pets or other objects before operating it safely.

Attention

When moving, please note:

- If you need to move it, it is recommended that two adults support the bottom of the armrest to lift the chair to move.
- Do not move the chair when it is unfolded to avoid damage to the chair parts.
- Do not directly push, pull and use other ways to move the chair to avoid scratching the floor and carpet or scuffing the footrest.
- Do not carry people or objects, and do not allow children to hand le it to avoid danger.

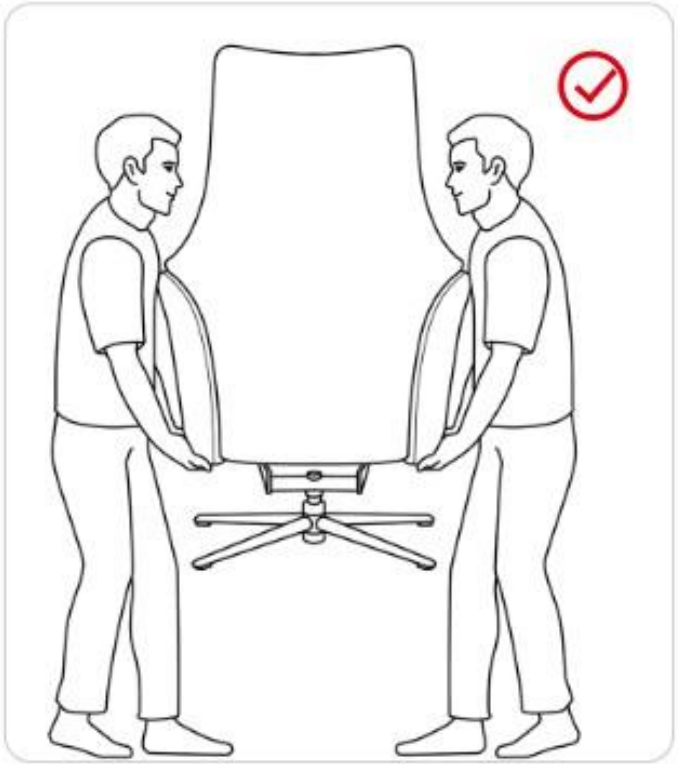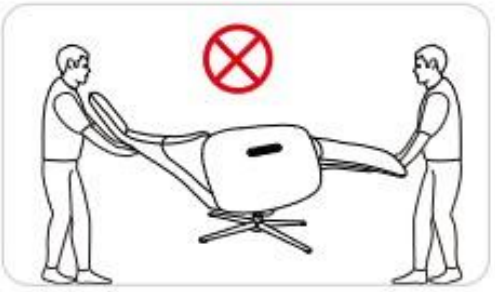Maintenance procedure:

· Daily maintenance procedure:

Use a clean soft cloth to clean fabric dust or stains frequently. The gap at the junction of the back, armrest and seat surface can be cleaned with a vacuum cleaner.

Attention

Special attention:

- Please place the chair in a place with a suitable temperature to avoid staining the fabric by moisture and dust, and avoid discoloration and aging of the fabric caused by strong light.
- Be sure to pull off the power wire 's power plug from the power outlet when performing maintenance, and be careful not to pull it out with wet hands to avoid electric shock and injury
- It is forbidden to use wet cloth, hard objects or chemicals such as acid and alkaline to touch the fabric, so as not to affect the surface quality and service cycle.
- Do not iron the fabric to avoid damage to the fabric.

· Regular maintenance procedure:

Generally, the maintenance shall be carried out once a quarter. After the chair surface is wiped with a clean soft cloth, the leather fabric can be evenly coated with special leather cleaner or polishing agent, and then cleaned with a clean soft cloth later. When it is not in use for a long time, please put on a protective cover to prevent dust on the surface. .

Attention

Special attention:

- Avoid children stabbing with sharp objects, and avoid children jumping and playing on the chair to avoid injury.
- Children should avoid scribbling on the chair with colored pens, markers or color powder and other pollutants. Once the pollution is found, it must be cleaned up in time.
- Children should not play with the power supply, pull the power wire and remove the battery pack alone.
- When sitting or lying on the sofa, it is not advisable to carry hard objects or metal products on the body, so as to avoid scratching the surface and making it difficult to repair, thus affecting the beauty.

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